

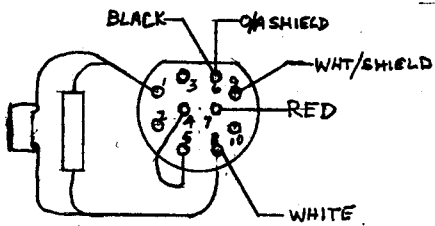



WIRING DIAGRAM

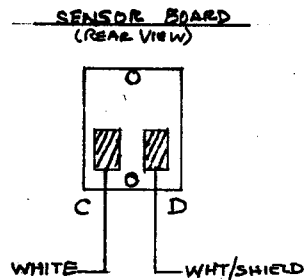

REAR VIEW

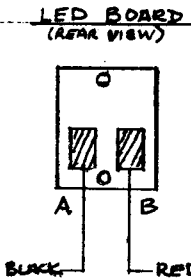

SCHEMATIC

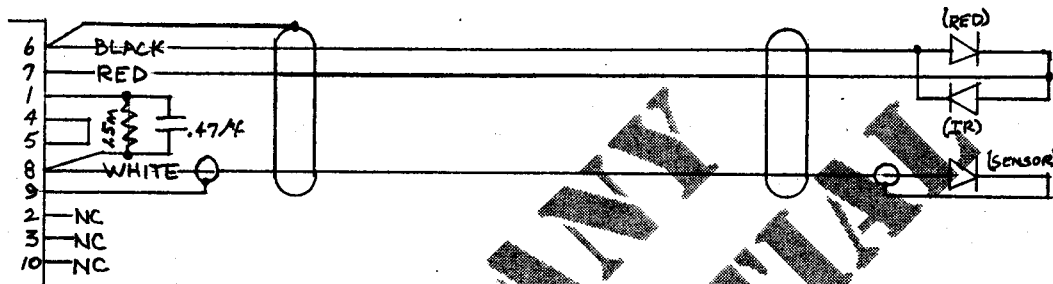

PARTS LIST

12' X2D05580
X3B01505
X3A04703
X3C11719
7 X3C20106

CABLE
RESISTOR 1/2W 1%
CAPACITOR 470P
CONNECTOR FEMALE, 10 PIN
CONNECTOR PINS, FEMALE

SPECIAL INSTRUCTIONS

1. WHEN REPLACING CABLE USE THE FOLLOWING WIRING DIAGRAM:

RED - PIN 6 TO A
YELLOW - 7 TO B
WHITE - 8 TO C
BLACK - NC
W/B SHIELD - 9 TO D
O/A SHIELD - NC



EPIC MEDICAL EQUIPMENT SERVICES, INC.

Dallas, Texas

SCALE: NA

APPROVED BY:

DRAWN BY WORLEY

DATE: 11-27-95

REVISÉ

DATEX SP02 CABLE - #878579-2 } OLD
SAS-F4 } VERSIONS

REPAIR STANDARDS

DRAWING NUMBER

DATEXØ1