


PH

MASSIMO EXTENSION

10.10.01.


~~CHANGING~~ cable original cable

Contacts FACE DOWN

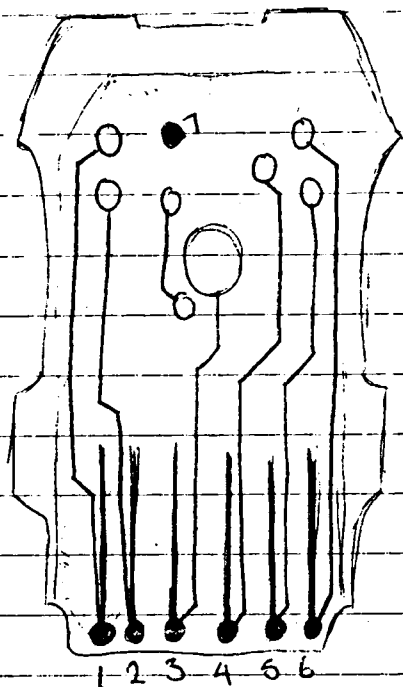

- | | |
|------------------------------|------------------|
| 1/ RED | 1/ YELLOW |
| 2/ RED ORANGE | 2/ RED |
| 3/ MAINSHIELD N/C | 3/ MS |
| 4/ INNERSHIELD | 4/ IS |
| 5/ WHITE | 5/ WHITE |
| 6/ GREEN | 6/ BLACK |
| 7/ MAINSHIELD | 7/ MS |

TOP VIEW

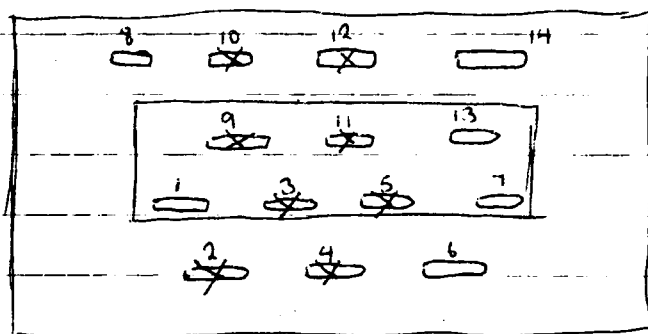

1/ RED

2/ N/C

3/ N/C

4/ N/C

5/ N/C

6/ Innershield

7/ GREEN

8/ ORANGE

9/ N/C

10/ N/C

11/ N/C

12/ N/C

13/ MAINSHIELD + CASING.

14/ WHITE

1/ YELLOW

2/

3/

4/

5/

6/

7/ BLACK

8/ RED

9/

10/

11/

12/


13/ MS + CASE

14/ WHITE

MASSIMO REPAIR Changing cable

Connector TOP VIEW

- 1/ Black
- 2/ white
- 3/ Inner shield
- 4/ outer shield
- 5/ ~~orange~~ Red
- 6/ ~~Red~~ Yellow

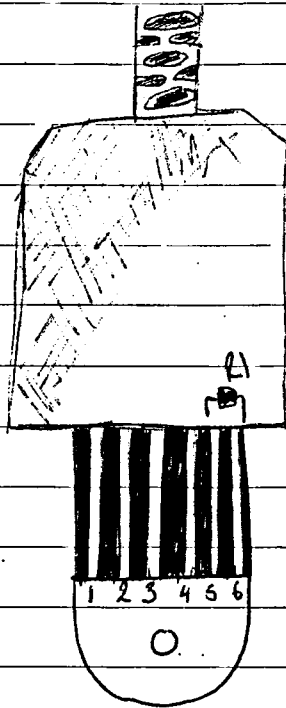


$R_1 = 12.72 \text{ k}\Omega$ in connector

LED F/V



- 1/ Red
- 2/ Yellow


SENSOR F/V

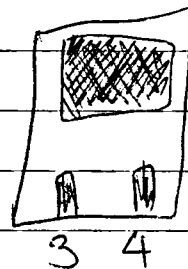

- 3/ White
- 4/ Black


Connector 2 wiring

Pin 1
Pin 2
Pin 3
Pin 4
Pin 5
Pin 6
Pin 7
Pin 8
Pin 9

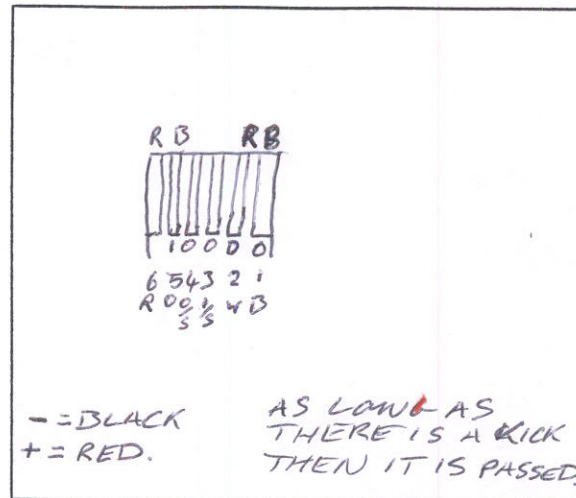

Test Instructions

Test using Test fixture curve tracer

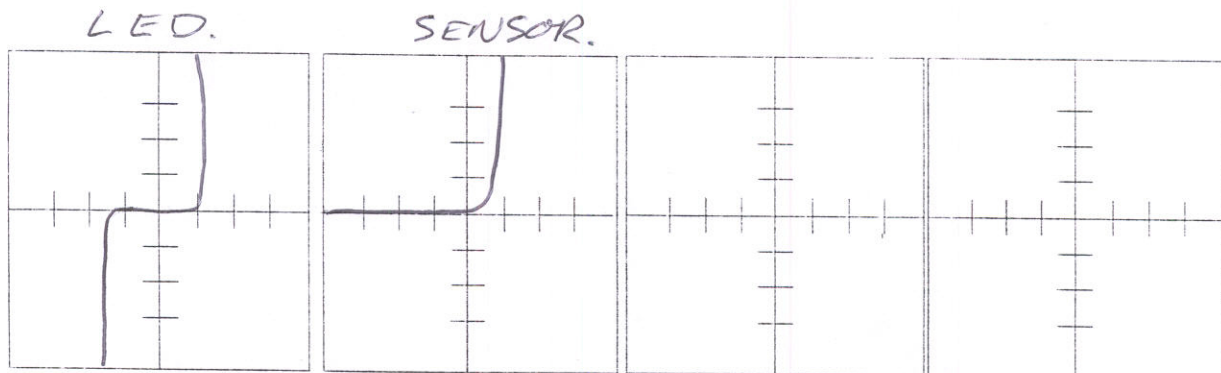

Drawn BY	
Date	
Checked By	
Date	
Revised By	
Revision Date	
Revision number	