

Pulse Oximeter Monitor: BCI 71000A2
Monitor S/N: 170276621

Finger Probe Model: P861RA
BioTek Simulator S/N: 125672

9731-2	7H00195	98	94	80	70	60	0	0	0	0	0
9731-2	7H00196	98	94	80	70	60	0	0	0	0	0
9731-3	7H00197	98	94	80	70	60	0	0	0	0	0
9731-3	7H00198	98	94	80	70	60	0	0	0	0	0
9731-3	7H00199	98	94	80	70	60	0	0	0	0	0
9731-3	7H00200	98	94	80	70	60	0	0	0	0	0
9731-3	7H00201	98	94	80	70	60	0	0	0	0	0
9731-3	7H00202	98	94	80	70	60	0	0	0	0	0
9731-3	7H00203	98	94	80	70	60	0	0	0	0	0
9731-3	7H00204	98	94	80	70	60	0	0	0	0	0
9731-3	7H00205	98	94	80	70	60	0	0	0	0	0
9731-3	7H00206	98	94	80	70	60	0	0	0	0	0
9731-3	7H00207	98	94	80	70	60	0	0	0	0	0
9731-3	7H00208	98	94	80	70	60	0	0	0	0	0
9731-3	7H00209	98	94	80	70	60	0	0	0	0	0
9731-3	7H00210	98	94	80	70	60	0	0	0	0	0
9731-3	7H00211	98	94	80	70	60	0	0	0	0	0
9731-3	7H00212	98	94	80	70	60	0	0	0	0	0
9731-3	7H00213	98	94	80	70	60	0	0	0	0	0
9731-3	7H00214	98	94	80	70	60	0	0	0	0	0
9731-3	7H00215	98	94	80	70	60	0	0	0	0	0
9731-3	7H00216	98	94	80	70	59	0	0	0	0	-1
9731-3	7H00217	98	94	80	70	60	0	0	0	0	0
9731-3	7H00218	98	94	80	70	60	0	0	0	0	0
9731-3	7H00219	98	94	80	70	60	0	0	0	0	0
9731-3	7H00220	98	94	80	70	59	0	0	0	0	-1
9731-3	7H00221	98	94	81	71	61	0	0	1	1	1
9731-3	7H00222	98	94	80	70	60	0	0	0	0	0
9731-3	7H00223	98	94	80	70	60	0	0	0	0	0
9731-3	7H00224	98	94	80	70	60	0	0	0	0	0
9731-3	7H00225	98	94	80	70	60	0	0	0	0	0
9731-3	7H00226	98	94	80	70	60	0	0	0	0	0
9731-3	7H00227	98	94	80	70	60	0	0	0	0	0
9731-3	7H00228	98	94	80	70	60	0	0	0	0	0
9731-3	7H00229	98	94	80	70	60	0	0	0	0	0
9731-3	7H00230	98	94	80	70	60	0	0	0	0	0
9731-3	7H00231	98	94	80	70	60	0	0	0	0	0
9731-3	7H00232	98	94	80	70	60	0	0	0	0	0
9731-3	7H00233	98	94	80	70	60	0	0	0	0	0
9731-3	7H00234	98	94	80	70	60	0	0	0	0	0
9731-3	7H00235	98	94	80	70	60	0	0	0	0	0
9731-3	7H00236	98	94	80	70	60	0	0	0	0	0
9731-3	7H00237	98	94	80	70	60	0	0	0	0	0
9731-3	7H00238	98	94	81	71	61	0	0	1	1	1
9731-3	7H00239	98	94	80	70	60	0	0	0	0	0
9731-3	7H00240	98	94	80	70	60	0	0	0	0	0
9731-3	7H00241	98	94	80	70	60	0	0	0	0	0
9731-3	7H00242	98	94	80	70	60	0	0	0	0	0
9731-3	7H00243	98	94	80	70	60	0	0	0	0	0
9731-3	7H00244	98	94	80	70	60	0	0	0	0	0
9731-3	7H00245	98	94	80	70	60	0	0	0	0	0

Tester: *Patricia Lee*

Date: 9/23/97

Q.A.: *JB*

Date: 9/23/97

A6.2 Test Results P861RA
Pulse Oximeter Monitor: BCI 71000A2
Monitor S/N: 170276621

Pulse Rate Testing
Finger Probe Model: P861RA
BioTek Simulator S/N: 125672

Lot Number	Serial Number	Test Results					Bias Error				
		250	150	95	80	50	250	150	95	80	50
9731-2	7H00147	250	150	95	80	50	0	0	0	0	0
9731-2	7H00148	250	150	95	80	50	0	0	0	0	0
9731-2	7H00149	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00150	250	150	95	80	50	0	0	0	0	0
9731-2	7H00151	250	150	95	80	50	0	0	0	0	0
9731-2	7H00152	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00153	250	150	95	80	50	0	0	0	0	0
9731-2	7H00154	250	150	95	80	50	0	0	0	0	0
9731-2	7H00155	250	150	95	80	50	0	0	0	0	0
9731-2	7H00156	250	150	95	80	50	0	0	0	0	0
9731-2	7H00157	250	150	95	80	50	0	0	0	0	0
9731-2	7H00158	250	150	95	80	50	0	0	0	0	0
9731-2	7H00159	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00160	250	150	95	80	50	0	0	0	0	0
9731-2	7H00161	250	150	95	80	50	0	0	0	0	0
9731-2	7H00162	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00163	250	150	95	80	50	0	0	0	0	0
9731-2	7H00164	250	150	95	80	50	0	0	0	0	0
9731-2	7H00165	250	150	95	80	50	0	0	0	0	0
9731-2	7H00166	250	150	95	80	50	0	0	0	0	0
9731-2	7H00167	250	150	95	80	50	0	0	0	0	0
9731-2	7H00168	250	150	95	80	50	0	0	0	0	0
9731-2	7H00169	250	150	95	80	50	0	0	0	0	0
9731-2	7H00170	250	150	95	80	50	0	0	0	0	0
9731-2	7H00171	250	150	95	80	50	0	0	0	0	0
9731-2	7H00172	250	150	95	80	50	0	0	0	0	0
9731-2	7H00173	250	150	95	80	50	0	0	0	0	0
9731-2	7H00174	250	150	95	80	50	0	0	0	0	0
9731-2	7H00175	250	150	95	80	50	0	0	0	0	0
9731-2	7H00176	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00177	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00178	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00179	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00180	250	150	95	80	50	0	0	0	0	0
9731-2	7H00181	250	150	95	80	50	0	0	0	0	0
9731-2	7H00182	250	150	95	80	50	0	0	0	0	0
9731-2	7H00183	250	150	95	80	50	0	0	0	0	0
9731-2	7H00184	250	150	95	80	50	0	0	0	0	0
9731-2	7H00185	250	150	95	80	50	0	0	0	0	0
9731-2	7H00186	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00187	250	150	95	80	50	0	0	0	0	0
9731-2	7H00188	249	150	95	80	50	-1	0	0	0	0
9731-2	7H00189	250	150	95	80	50	0	0	0	0	0
9731-2	7H00190	250	150	95	80	50	0	0	0	0	0
9731-2	7H00191	250	150	95	80	50	0	0	0	0	0
9731-2	7H00192	250	150	95	80	50	0	0	0	0	0
9731-2	7H00193	250	150	95	80	50	0	0	0	0	0
9731-2	7H00194	250	150	95	80	50	0	0	0	0	0

Pulse Oximeter Monitor: BCI 71000A2
Monitor S/N: 170276621

Finger Probe Model: P861RA
BioTek Simulator S/N: 125672

9731-2	7H00195	250	150	95	80	50	0	0	0	0	0
9731-2	7H00196	250	150	95	80	50	0	0	0	0	0
9731-3	7H00197	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00198	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00199	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00200	250	150	95	80	50	0	0	0	0	0
9731-3	7H00201	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00202	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00203	250	150	95	80	50	0	0	0	0	0
9731-3	7H00204	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00205	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00206	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00207	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00208	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00209	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00210	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00211	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00212	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00213	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00214	250	150	95	80	50	0	0	0	0	0
9731-3	7H00215	250	150	95	80	50	0	0	0	0	0
9731-3	7H00216	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00217	250	150	95	80	50	0	0	0	0	0
9731-3	7H00218	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00219	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00220	250	150	95	80	50	0	0	0	0	0
9731-3	7H00221	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00222	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00223	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00224	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00225	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00226	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00227	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00228	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00229	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00230	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00231	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00232	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00233	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00234	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00235	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00236	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00237	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00238	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00239	250	150	95	80	50	0	0	0	0	0
9731-3	7H00240	250	150	95	80	50	0	0	0	0	0
9731-3	7H00241	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00242	250	150	95	80	50	0	0	0	0	0
9731-3	7H00243	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00244	249	150	95	80	50	-1	0	0	0	0
9731-3	7H00245	249	150	95	80	50	-1	0	0	0	0

Tester: *[Signature]*

Date: 9/23/97

Q.A.: *[Signature]*

Date: 9/23/97

A6.3 P862RA Test Results

Oximeter Monitor: CSI 504

Monitor S/N: 291108486

SpO₂ Testing

Finger Probe Model: P862RA

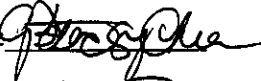
BioTek Simulator S/N: 125672

Lot Number	Serial Number	Test Results					Bias Error				
		98	94	80	70	60	98	94	80	70	60
9731-2	7H00880	99	95	82	71	62	1	1	2	1	2
9731-2	7H00881	98	95	82	70	62	0	1	2	0	2
9731-2	7H00882	99	96	82	71	62	1	2	2	1	2
9731-2	7H00883	98	96	81	71	61	0	2	1	1	1
9731-2	7H00884	99	96	82	71	62	1	2	2	1	2
9731-2	7H00885	98	96	80	70	62	0	2	0	0	2
9731-2	7H00886	98	95	80	70	60	0	1	0	0	0
9731-2	7H00887	98	96	82	71	61	0	2	2	1	1
9731-2	7H00888	99	96	82	71	62	1	2	2	1	2
9731-2	7H00889	98	96	81	71	61	0	2	1	1	1
9731-2	7H00890	98	96	82	71	61	0	2	2	1	1
9731-2	7H00891	98	96	82	71	61	0	2	2	1	1
9731-2	7H00893	99	96	80	70	60	1	2	0	0	0
9731-2	7H00894	99	96	82	70	60	1	2	2	0	0
9731-2	7H00895	99	96	81	71	62	1	2	1	1	2
9731-2	7H00897	98	95	82	70	62	0	1	2	0	2
9731-2	7H00898	98	96	81	71	60	0	2	1	1	0
9731-2	7H00899	98	96	82	71	61	0	2	2	1	1
9731-2	7H00900	98	95	80	71	60	0	1	0	1	0
9731-2	7H00901	99	96	82	71	62	1	2	2	1	2
9731-2	7H00902	99	96	80	70	60	1	2	0	0	0
9731-2	7H00903	99	96	82	71	62	1	2	2	1	2
9731-2	7H00904	98	95	81	71	60	0	1	1	1	0
9731-2	7H00905	99	95	81	70	60	1	1	1	0	0
9731-2	7H00906	98	96	82	71	60	0	2	2	1	0
9731-2	7H00907	98	96	81	71	61	0	2	1	1	1
9731-2	7H00908	98	96	81	70	60	0	2	1	0	0
9731-2	7H00909	98	95	80	70	60	0	1	0	0	0
9731-2	7H00910	98	95	82	71	61	0	1	2	1	1
9731-2	7H00911	98	95	82	71	61	0	1	2	1	1
9731-2	7H00912	99	95	82	71	60	1	1	2	1	0
9731-2	7H00913	98	96	82	74	62	0	2	2	4	2
9731-2	7H00914	99	96	82	70	61	1	2	2	0	1
9731-2	7H00915	98	95	82	71	61	0	1	2	1	1
9731-2	7H00916	98	96	82	70	60	0	2	2	0	0
9731-2	7H00917	99	96	82	71	62	1	2	2	1	2

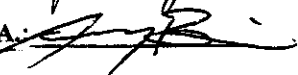
Oximeter Monitor: CSI 504
Monitor S/N: 291108486

Finger Probe Model: P862RA
BioTek Simulator S/N: 125672

9731-1	7I01432	98	96	81	70	61	0	2	1	0	1
9731-1	7I01433	98	96	82	71	62	0	2	2	1	2
9731-1	7I01434	98	96	82	72	61	0	2	2	2	1
9731-1	7I01435	98	96	82	72	61	0	2	2	2	1
9731-1	7I01436	98	96	80	69	59	0	2	0	-1	-1
9731-1	7I01437	98	96	82	72	62	0	2	2	2	2
9731-1	7I01438	98	94	80	70	60	0	0	0	0	0
9731-1	7I01439	98	96	82	71	61	0	2	2	1	1
9731-1	7I01440	98	96	82	72	61	0	2	2	2	0
9731-1	7I01441	99	96	81	70	60	1	2	1	0	-1
9731-1	7I01442	98	95	80	69	59	0	1	0	-1	2
9731-1	7I01443	99	96	83	72	62	1	2	3	2	0
9731-1	7I01444	98	95	81	70	60	0	1	1	0	-1
9731-1	7I01445	99	96	81	71	59	1	2	1	1	2
9731-1	7I01446	98	96	82	72	62	0	2	2	2	2
9731-1	7I01447	98	96	82	71	62	0	2	2	1	2
9731-1	7I01448	99	96	80	70	62	1	2	0	0	1
9731-1	7I01449	98	96	80	70	61	0	2	0	0	-1
9731-1	7I01450	98	96	80	69	59	0	2	0	-1	0
9731-1	7I01451	99	96	80	70	60	1	2	0	0	1
9731-1	7I01452	98	95	82	72	61	0	1	2	2	2
9731-1	7I01453	98	96	82	71	62	0	2	2	1	0
9731-1	7I01454	98	95	80	69	60	0	1	0	-1	0
9731-1	7I01455	98	96	80	72	60	0	2	0	2	0
9731-1	7I01456	98	96	80	70	60	0	2	0	0	0
9731-1	7I01457	96	95	80	70	60	-2	1	0	0	1
9731-1	7I01458	96	94	82	72	61	-2	0	2	2	1
9731-1	7I01459	95	94	82	71	61	-3	0	2	1	1
9731-1	7I01460	95	94	80	71	61	-3	0	0	1	2
9731-1	7I01461	95	94	82	72	62	-3	0	2	2	2
9731-1	7I01462	95	94	80	70	60	-3	0	0	0	0
9731-1	7I01463	96	94	81	71	62	-2	0	1	1	2
9731-1	7I01464	95	94	82	72	61	-3	0	2	2	1
9731-1	7I01465	95	94	81	70	60	-3	0	1	0	0
9731-1	7I01466	96	94	82	71	61	-2	0	2	1	1
9731-1	7I01467	96	94	83	72	61	-2	0	3	2	1
9731-1	7I01468	96	94	82	72	62	-2	0	2	2	2

Tester: 

Date: 9/24/97

Q.A.: 

Date: 9/24/97

A6.3 P862RA Test Results

Oximeter Monitor: CSI 504

Monitor S/N: 291108486

Pulse Rate Testing

Finger Probe Model: P862RA

BioTek Simulator S/N: 125672

Lot Number	Serial Number	Test Results					Bias Error				
		250	150	95	80	50	250	150	95	80	50
9731-2	7H00880	244	149	94	79	49	-6	-1	-1	-1	-1
9731-2	7H00881	244	149	94	79	49	-6	-1	-1	-1	-1
9731-2	7H00882	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00883	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00884	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00885	244	149	94	79	49	-6	-1	-1	-1	-1
9731-2	7H00886	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00887	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00888	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00889	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00890	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00891	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00893	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00894	243	149	94	80	49	-7	-1	-1	0	-1
9731-2	7H00895	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00897	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00898	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00899	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00900	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00901	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00902	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00903	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00904	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00905	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00906	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00907	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00908	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00909	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00910	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00911	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00912	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00913	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00914	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00915	243	149	94	79	49	-7	-1	-1	-1	-1
9731-2	7H00916	242	149	94	79	49	-8	-1	-1	-1	-1
9731-2	7H00917	242	149	94	79	49	-8	-1	-1	-1	-1

Oximeter Monitor: CSI 504
Monitor S/N: 291108486

Finger Probe Model: P862RA
BioTek Simulator S/N: 125672

9731-1	7I01432	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01433	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01434	241	149	94	79	49	-9	-1	-1	-1	-1
9731-1	7I01435	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01436	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01437	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01438	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01439	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01440	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01441	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01442	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01443	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01444	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01445	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01446	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01447	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01448	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01449	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01450	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01451	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01452	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01453	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01454	244	149	94	79	49	-6	-1	-1	-1	-1
9731-1	7I01455	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01456	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01457	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01458	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01459	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01460	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01461	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01462	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01463	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01464	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01465	242	149	94	79	49	-8	-1	-1	-1	-1
9731-1	7I01466	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01467	243	149	94	79	49	-7	-1	-1	-1	-1
9731-1	7I01468	252	149	94	79	49	2	-1	-1	-1	-1

Tester: [Signature]

Date: 9/26/97

Q.A.: [Signature]

Date: 9/26/97

A6.2 Test Results P861RA
Pulse Oximeter Monitor: BCI 71000A2
Monitor S/N: 170276621

SpO₂ Testing
Finger Probe Model: P861RA
BioTek Simulator S/N: 125672

Lot Number	Serial Number	Test Results					Bias Error				
		98	94	80	70	60	98	94	80	70	60
9731-2	7H00147	98	94	80	70	60	0	0	0	0	0
9731-2	7H00148	98	94	80	70	60	0	0	0	0	0
9731-2	7H00149	98	94	80	70	60	0	0	0	0	0
9731-2	7H00150	98	94	80	70	60	0	0	0	0	0
9731-2	7H00151	98	94	80	70	60	0	0	0	0	0
9731-2	7H00152	98	94	80	70	60	0	0	0	0	0
9731-2	7H00153	98	94	80	70	60	0	0	0	0	0
9731-2	7H00154	98	94	80	70	60	0	0	0	0	0
9731-2	7H00155	98	94	80	70	60	0	0	0	0	0
9731-2	7H00156	98	94	80	70	60	0	0	0	0	0
9731-2	7H00157	98	94	80	70	60	0	0	0	0	0
9731-2	7H00158	98	94	80	70	60	0	0	0	0	0
9731-2	7H00159	98	94	80	70	60	0	0	0	0	0
9731-2	7H00160	98	94	80	70	60	0	0	0	0	0
9731-2	7H00161	98	94	80	70	60	0	0	0	0	0
9731-2	7H00162	98	94	80	70	60	0	0	0	0	0
9731-2	7H00163	98	94	80	70	60	0	0	0	0	0
9731-2	7H00164	98	94	80	70	60	0	0	0	0	0
9731-2	7H00165	98	94	80	70	60	0	0	0	0	0
9731-2	7H00166	98	94	80	70	60	0	0	0	0	0
9731-2	7H00167	98	94	80	70	60	0	0	0	0	0
9731-2	7H00168	98	94	80	70	60	0	0	0	0	0
9731-2	7H00169	98	94	80	70	60	0	0	0	0	0
9731-2	7H00170	98	94	80	70	60	0	0	0	0	0
9731-2	7H00171	98	94	80	70	60	0	0	0	0	0
9731-2	7H00172	98	94	80	70	60	0	0	0	0	0
9731-2	7H00173	98	94	80	70	60	0	0	0	0	0
9731-2	7H00174	98	94	80	70	60	0	0	0	0	0
9731-2	7H00175	98	94	80	70	60	0	0	0	0	0
9731-2	7H00176	98	94	80	70	60	0	0	0	0	0
9731-2	7H00177	98	94	80	70	60	0	0	0	0	0
9731-2	7H00178	98	94	80	70	60	0	0	0	0	0
9731-2	7H00179	98	94	80	70	60	0	0	0	0	0
9731-2	7H00180	98	94	80	70	60	0	0	0	0	0
9731-2	7H00181	98	94	80	70	60	0	0	0	0	0
9731-2	7H00182	98	94	80	70	60	0	0	0	0	0
9731-2	7H00183	98	94	80	70	60	0	0	0	0	0
9731-2	7H00184	98	94	80	70	60	0	0	0	0	0
9731-2	7H00185	98	94	80	70	60	0	0	0	0	0
9731-2	7H00186	98	94	80	70	60	0	0	0	0	0
9731-2	7H00187	98	94	80	70	60	0	0	0	0	0
9731-2	7H00188	98	94	80	70	60	0	0	0	0	0
9731-2	7H00189	98	94	80	70	60	0	0	0	0	0
9731-2	7H00190	98	94	80	70	60	0	0	0	0	0
9731-2	7H00191	98	94	80	70	60	0	0	0	0	0
9731-2	7H00192	98	94	80	70	60	0	0	0	0	0
9731-2	7H00193	98	94	80	71	60	0	0	0	1	0
9731-2	7H00194	98	94	80	70	60	0	0	0	0	0