


P887 RA (Space Labs)


MS

Connector

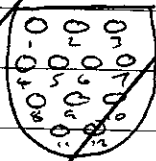

- 1) Resistor
- 2) Black
- 3) Resistor + Shields white
- 4) -
- 5) -
- 6) -
- 7) Resistor + Shields
- 8) Orange
- 9) Red
- 10) -
- 11) -
- 12) Yellow


Nellcor Extension Cable To Drager. (Brown/white Lemo).

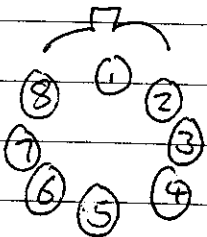

B

- 1) White
- 2) Yellow
- 3) Blue
- 4) Red
- 5) Black
- 6) - n/c
- 7) - n/c
- 8) Inner shield

R/V

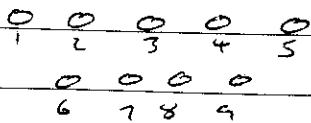

- 1) Blue
- 2) Red
- 3) Black
- 4) - n/c
- 5) White
- 6) Yellow
- 7) Main Shield
- 8) - n/c
- 9) Inner shield.


R/V

Date	24/03/97	Type	Nellcor Original
Nellcor.cmx			
Sw1 Sw2			
Drawn By	D.Lamb	Signed	

rell cov

Rear view

Rie-camp

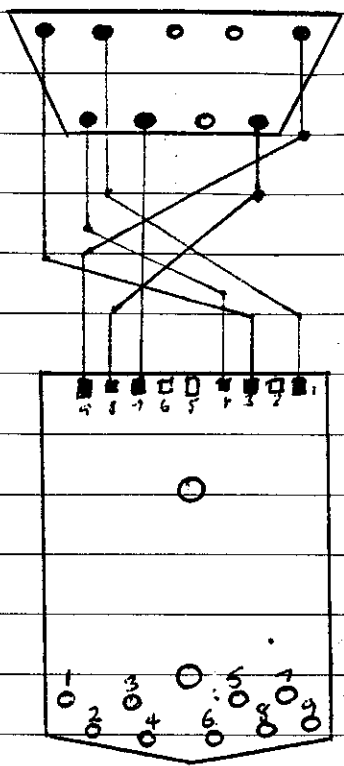

- 1 Blue
- 2 Black
- 3 Red
- 4 Brown
- 5 Green
- 6 Orange
- 7 Yellow
- 8 Purple
- 9 Grey

Date		Type	
------	--	------	--

NEELCOR EXT. CABLE

9 PIN DB FEMALE


NEELCOR PLUS
FRONT VIEW

I.R.	L.E.D.	Sensor
---	---	---

Drawn By:	Signed
-----------	--------

BCI - Nellcor extension


male end - RV

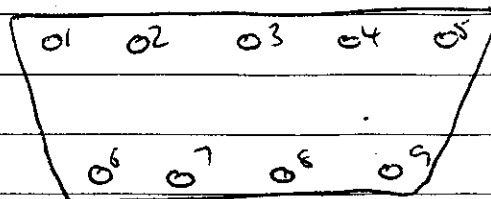

- 1 Black 6 White
- 2 n/c 7 n/c
- 3 Yellow 8 OS
- 4 Red 9 Light + IS
- 5 Resistor

THIS DOES NOT WORK

female end - RV

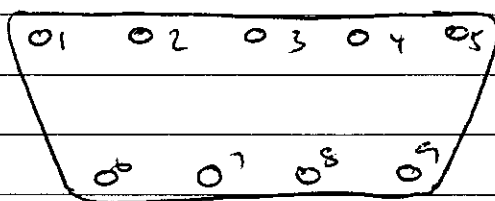


- 1 n/c 6 IS
- 2 Red 7 OS
- 3 Yellow 8 n/c
- 4 n/c 9 White
- 5 Black

BCI - Nellcor extension cable

USING VIAMED STOCK PSS8RA CABLES

P/N 0018587


Female end - RV

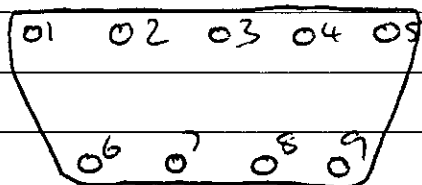

- 1 n/c
- 2 Yellow
- 3 Red
- 4 n/c
- 5 White
- 6 ~~Black~~ NC
- 7 Outer shield
- 8 n/c
- 9 Black

This works

PSS9RA with male conn. - PRELIM try! for evaluation by customer

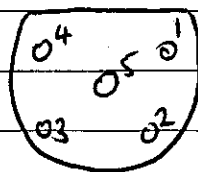

- 1 Yellow
- 2 Red
- 3 Black
- 4 White
- 5 Shields