

Disinfectant Test Report:

bluepoint MEDICAL

Confidential - Distribution only by written permission of bluepoint medical GmbH & Co. KG

Devices

Sensor Tested:	SoftCap	Monitor Tested:	Masimo Radical 7
Sensor REF:	6020132302	Monitor SN:	119115
Sensor SN:	120905-004		
Sensor Tested:	SoftWrap		
Sensor REF:	6020132310		
Sensor SN:	130902-005		
Sensor Tested:	SoftFlap		
Sensor REF:	6020132303		
Sensor SN:	130901-002		

Method

Simulator:	Index 2 (Bio-Tek)	R-Curve:	Masimo
Serial number:	158657	SpO2:	84%
Calibrated until:	August 2014		
Heart Rate:	70 bpm	Test performed:	R. Jürgensen
Disinfectant:	Isopropanol 70%, Meliseptol (B.Braun), Teccare control (Talley Environmental Care)		
Method:	Sensors are cleaned 20 times with disinfectant and a soft cloth		
Steady-State:	Measurement was taken after value was stable (5 - 10 sec)		

Test Results

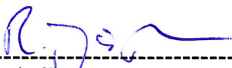
Isopropanol 70%	SoftCap	SoftWrap	SoftFlap
Visual test	<i>passed</i>	<i>passed</i>	<i>passed</i>
Functional test	<i>passed</i>	<i>passed</i>	<i>passed</i>

Meliseptol	SoftCap	SoftWrap	SoftFlap
Visual test	<i>passed</i>	<i>passed</i>	<i>passed</i>
Functional test	<i>passed</i>	<i>passed</i>	<i>passed</i>

Teccare control	SoftCap	SoftWrap	SoftFlap
Visual test	<i>passed</i>	<i>passed</i>	<i>passed</i>
Functional test	<i>passed</i>	<i>passed</i>	<i>passed</i>

Conclusion

The tested sensor meets the stated accuracy specification of +/- 3bpm at 70 bpm.


Rayk Jürgensen

Test Date(s): 18.12.2013

Report Date: 18.12.2013