


bluepoint MEDICAL	Project Disposable SPO2 sensors	Revision 0
Development Document	Labeling Nellcor compatible type 1-AP-VM (Viamed)	Page 1 of 3

Design Output

File No.: LO-01-10-0001 Rev.0
Date: 2015-04-27

Pouch insert label:

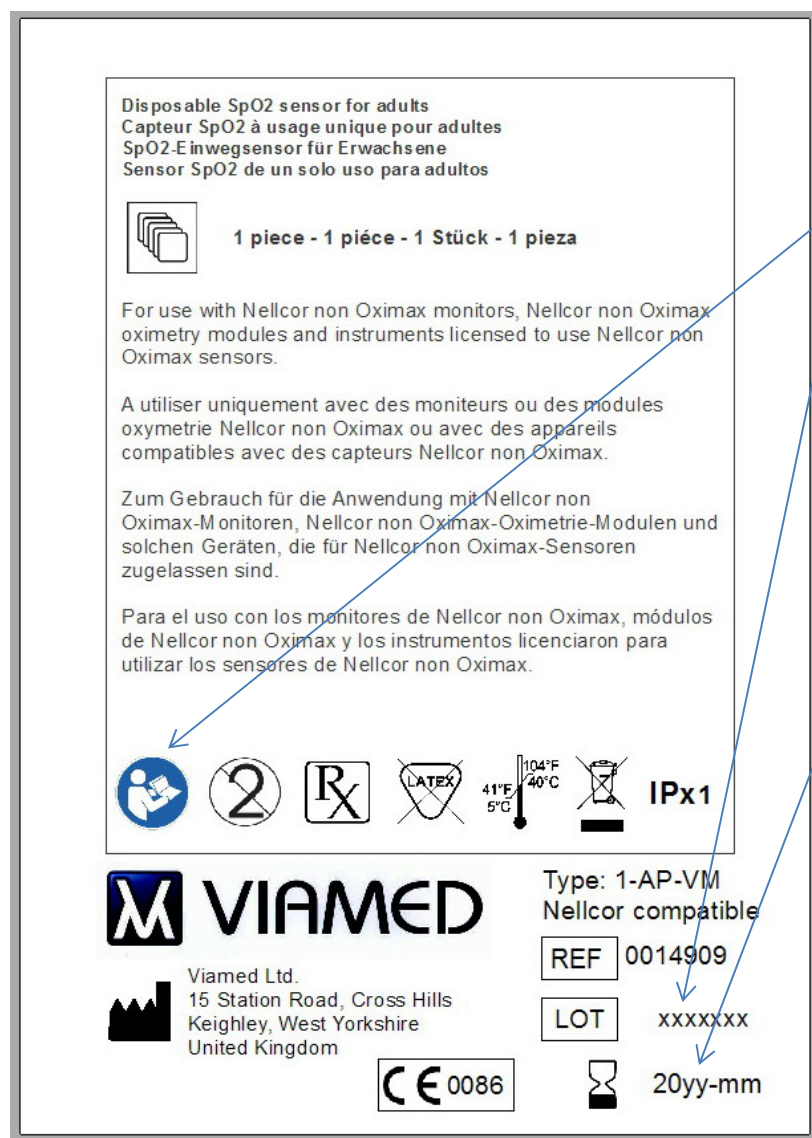

Requirements:

- 1) Insert material: White paper carton
- 2) Insert dimensions (L x W): 90 x 130 [mm]
- 3) Print: Black/White
- 4) Symbol ISO 7010-M002 blue print or blue sticker
- 5) LOT numbering:

X X X X X X X

→ Continuous numbers (001, 002, ...)

→ Sensor type I (I= 1-AP-VM)

→ Month of manufacturing (A= JAN, B= FEB, C= ...)


→ Year of manufacturing (A= 2015, B= 2016, C= ...)

→ Manufacturer 1

LOT Example: 1ADI123
- 6) Year and Month of applicability
Example: 2018-04

bluepoint MEDICAL	Project Disposable SPO2 sensors	Revision 0
Development Document	Labeling Nellcor compatible type 1-AP-VM (Viamed)	Page 2 of 3

Packaging Box top label:

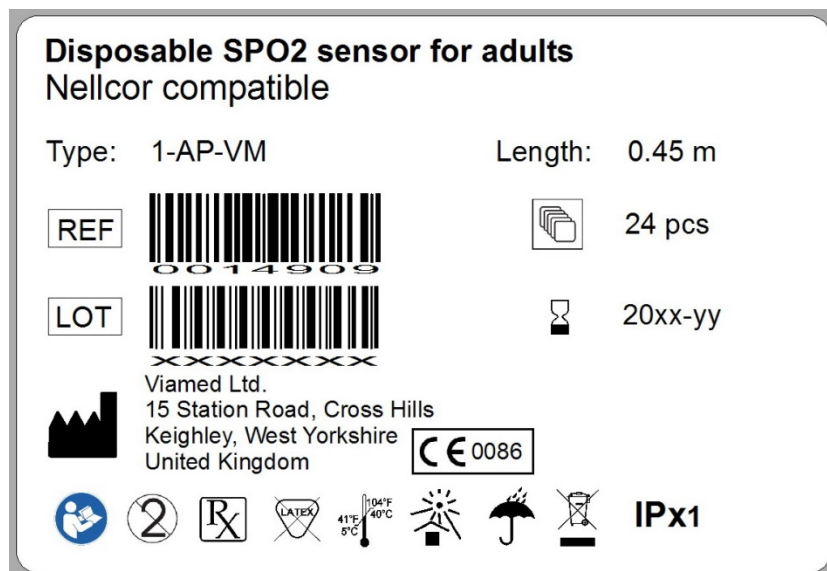

Requirements:

- 1) Label material: White semigloss paper
- 2) Label print color: Black/White
- 3) Symbol ISO 7010-M002 blue print or blue sticker
- 4) LOT numbering referred to insert requirements
- 5) Label dimensions (L x W): 100 x 80 [mm]
- 6) Label positioning: Centered on carton box top side

Packaging Box side label:



Requirements:

- 1) Label material: White semigloss paper
- 2) Label print color: Black/White
- 3) Label dimensions (L x W): 32 x 19 [mm]
- 4) Label positioning: Centered on carton box left side

	Project Disposable SPO2 sensors	Revision 0
Development Document	Labeling Nellcor compatible type 1-AP-VM (Viamed)	Page 3 of 3

Sensor Compatibility hang tag:

Nellcor compatible	For use with non OXIMAX technology monitors or modules
--------------------	--

Requirements:

- | | |
|------------------------------|--|
| 1) Label material: | White polyester |
| 2) Label print color: | Black/White |
| 3) Label dimensions (L x W): | 90 x 15 [mm] |
| 4) Label positioning: | 150mm from connector outer edge |
| 5) Readability: | Connector at the bottom, sensor on top |