



2305 South 1070 West
West Valley, UT 84119

DOCUMENT NUMBER: RP34P03-001-LCS

DOCUMENT TITLE: Label File, 0-15 LPM Flowmeter Green 50 PSI

DOCUMENT NOTES:

Document Information

Revision: D

Vault: Label-rel

Doc Type: Labeling

Status: Release

Date Information

Effective Date: 09 Jan 2024

Next Review Date: 09 Jan 2026

Release Date: 09 Jan 2024

Expiration Date:

Control Information

Author: DLARSEN

Previous Number:

Owner: DLARSEN

Change Number: DC-0372-23

All dates and times are in US/Mountain

LABEL CONTROL SHEET

Revision # D MAXTEC PART NUMBER: RP34P03-001 DESCRIPTION: 0-15 LPM FLOWMETER GREEN 50 PSI Num. of Labels: 4

GTIN NUMBER SHOULD READ: 00817770022664

LABEL USING WORK INSTRUCTIONS: WI-0037

Use FRM-0292

LABEL STOCK NUM: RP35P13-001
 LABEL DESCRIPTION: 4 X 3.5 Cardstock Label Stock
 FILE PATH: Z:\LABEL FILE\ FLOWMETER/ RP34P03-001

TO PRINT LABEL CHANGE LOT NUMBER USING RECEIVING LOT NUMBER, CHANGE MANUFACTURING DATE USING ORIGINAL LABEL DATE AND YYYY-MM-DD SEQUENCE CHANGE BARCODE LOT NUMBER BY GOING INTO PROPERTIES CHANGE THE DIGITS AFTER THE (10) TO THE LOT NUMBER, CHANGE THE DIGITS AFTER THE (11) AND BEFORE THE (10) USING DATE FROM ORIGINAL LABEL AND DATE CODE YYMMDD AND USE THERMAL TRANSFER RIBBON

FLOWMETER 0-15 LPM, GREEN 50 PSI		CONTENTS: 1 FLOWMETER 1 MANUAL	
REF RP34P03-001		EMERGO EUROPE Westervoortsedijk 60, 6827 AT Arnhem The Netherlands	
LOT 120121 2022-02-14			
1 Each	MD	Emergo Consulting (UK) Limited C/O CR 360 - UL International, Compass House, Vision Park Histon, Cambridge CB24 9BZ, United Kingdom	
140°F (60°C) 95 R _x ONLY	UK REP		
0123	Precision Medical, Inc. Northampton, PA 18067	Manufactured For: maxtec 2305 South 1070 West Salt Lake City, Ut 84119 USA 801-266-5300 1-800-748-5355 Fax: 801-270-5590	
MADE IN USA			
TP-0034 (02)		UDI	
(01)00817770022664(11)220214(10)120121			

MAXTEC

OPERATIONAL PROCEDURE

LABEL STOCK NUM: RP35P66

LABEL DESCRIPTION: 1 X .5 WHITE LABEL STOCK

FILE PATH: Z:\LABEL FILE FLOWMETER RP34P03-001 UDI

TO PRINT LABEL CHANGE MANUFACTURING DATE USING ORIGINAL LABEL DATE AND YYYY-MM-DD SEQUENCE CHANGE BARCODE LOT NUMBER BY GOING INTO DATA SOURCE PROPERTIES CHANGE THE DIGITS AFTER THE (10) TO THE LOT NUMBER, CHANGE THE DIGITS AFTER THE (11) AND BEFORE THE (10) USING ORIGINAL LABEL DATE AND DATE CODE YYMMDD AND USE THERMAL TRANSFER RIBBON




LABEL STOCK NUM: RP35P70

LABEL DESCRIPTION: .800 X .25 LABEL STOCK

FILE PATH: Z:\LABEL FILE FLOWMETER RP34P03-001 UNIT

TO PRINT LABEL CHANGE LOT USING RECEIVING LOT NUMBER AND USE THERMAL TRANSFER RIBBON

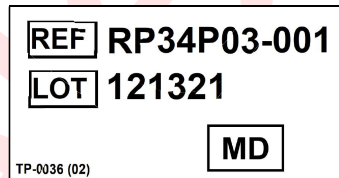

LABEL ON PRODUCT SHOULD BE:

