


REMOVE ALL SHARP EDGES

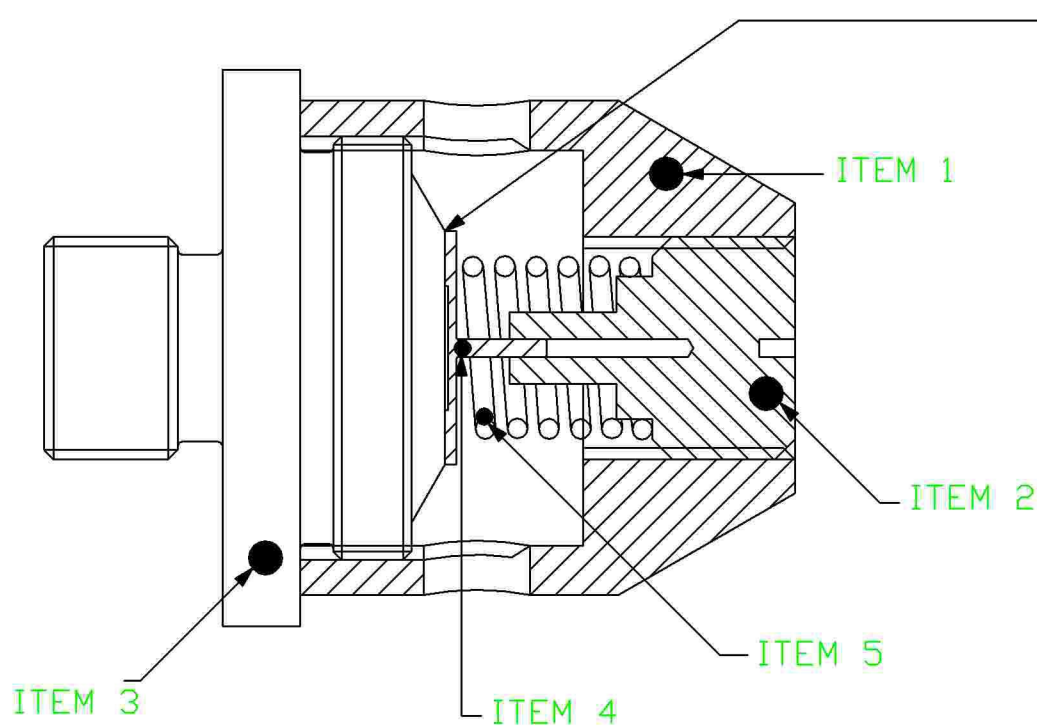

Disc - Diamond Lapped to
Optical Flat

Knife Edge - Diamond Lapped to
Optical Flat

Maximum width of seat - 0.008"

Set to 30cm H₂O at 5L/min Flow

Apply Kilopoise 0001G to Spindle

PARTS LIST

Item 1 - Valve Body	0330224
Item 2 - Adjustment Screw	0330228
Item 3 - Valve Seating	0330225
Item 4 - Valve Seat	0330227
Item 5 - Spring	0330222

Material : See Parts List					Title Precision Valve Assembly		 Dim in mm Tol ± 0.2		VIAMED Ltd. 15 Station Rd Cross Hills, Keighley West Yorkshire BD20 7DT	
A	13/02/04	1	K.J.Rush		Scale	Not To Scale	Drawn	K. Rush		
REV	Date	N°	Drawn	Approv	Part No	0330210	Date	22/01/04	Dwg No.	TT0100.05