


capnography


VersaStream[®]

Luer Lock Gas Sampling Lines

Innovative Solutions for Sidestream Gas Monitoring for use with
Capnographs and Multigas Monitors that incorporate a Luer Lock sampling port

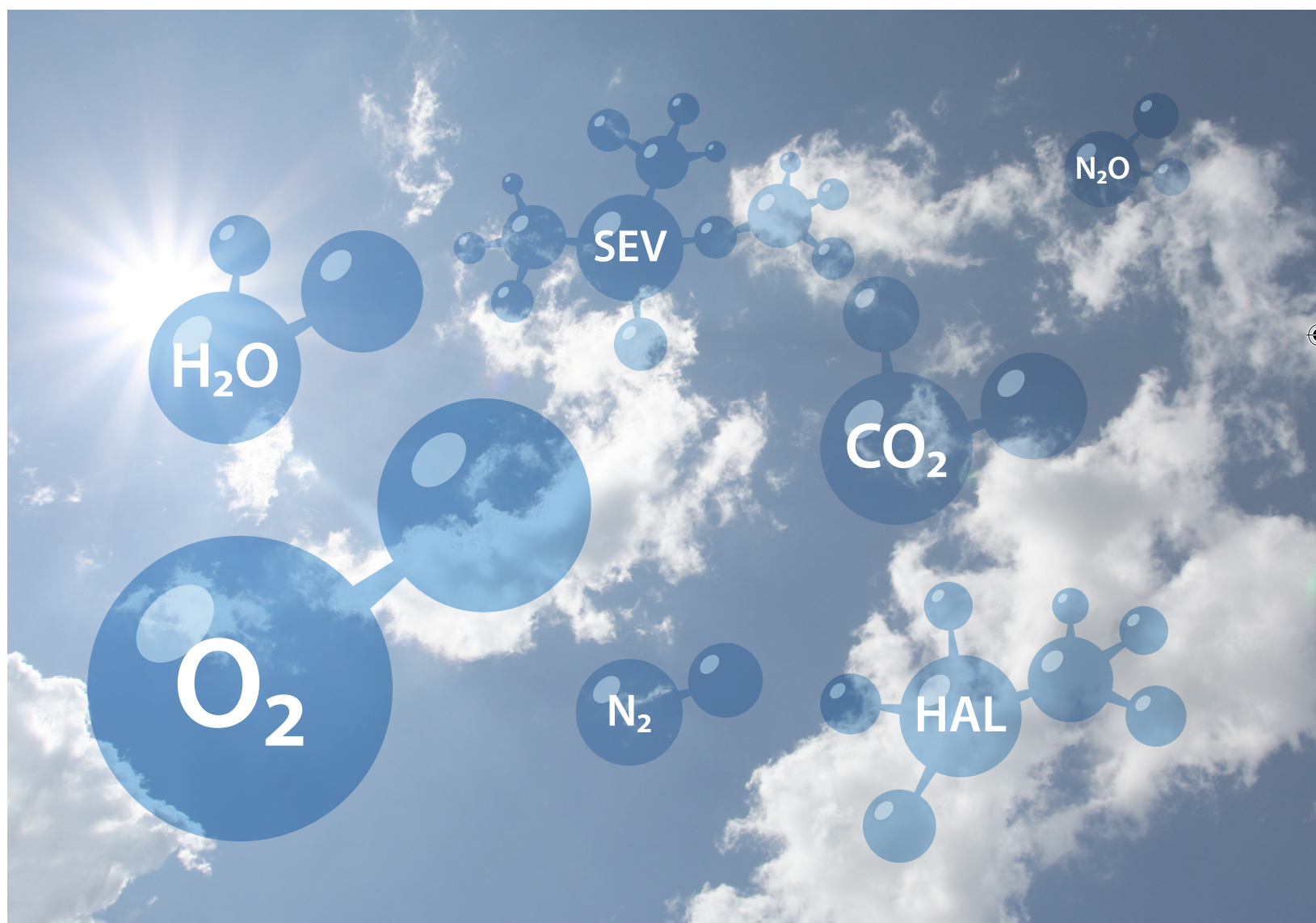

VersaStream® gas sampling lines

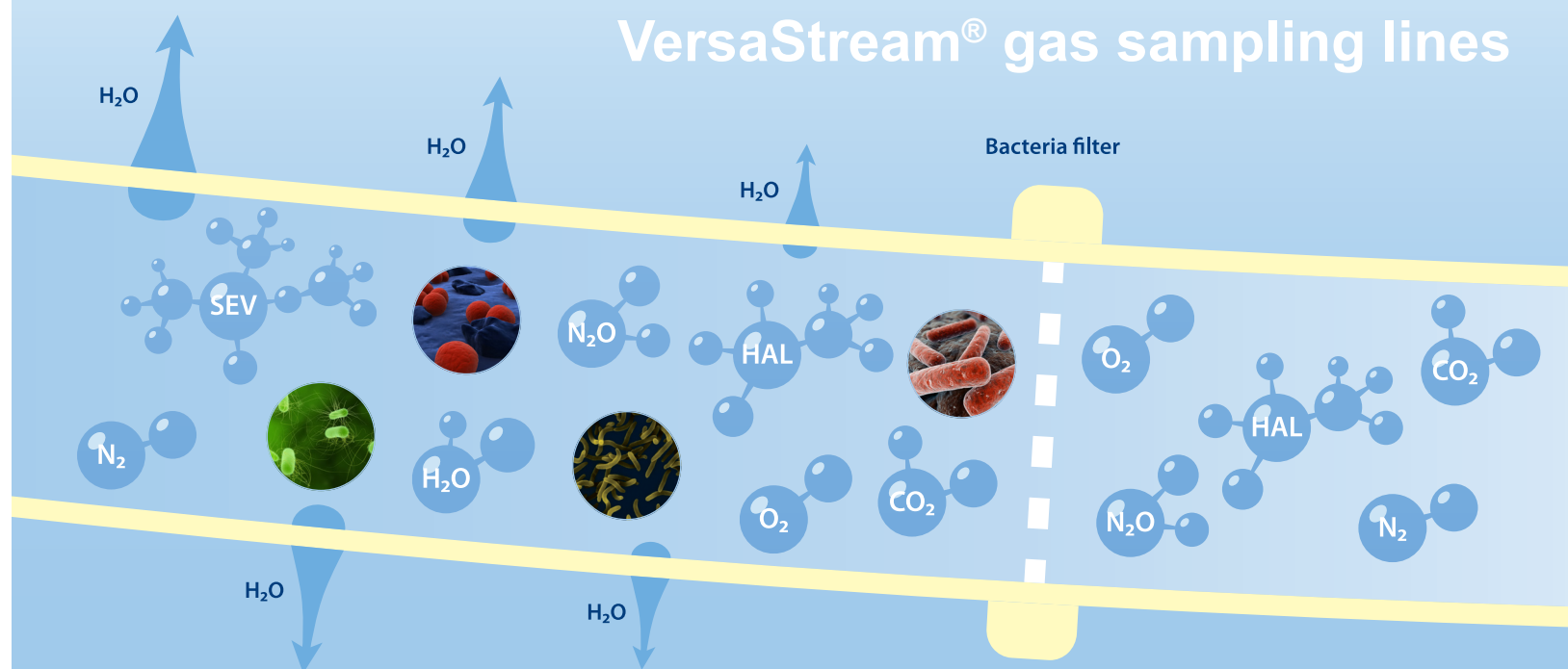

The VersaStream® gas sampling lines combine reliability, flexibility and cost efficiencies

Why VersaStream® gas sampling lines?


- Cost-effective alternative to Nafion® tubes.
- The gas-selective non-porous material facilitates fast and complete removal of condensation from the sampling line.
- The integrated bacterial filter prevents contamination.
- Regular water trap exchange (if present) is not required.
- Reliable measured gas values, due to negligible dead space volume of the filter - even in cases of high respiration rates.
- For use with many common types of gas exchange monitors, either low-flow or high-flow.
- Integrated bacteria filter and water removal reduces the number of patient circuit parts and speeds up patient monitoring preparation time.

Innovative material technologies

Latest material technologies permit active and complete removal of moisture and water – even in cases of low temperatures and high humidity.

No need to use expensive Nafion® tubes, which are not resistant to bending.


VersaStream® gas sampling lines are very flexible and are suitable for short-term and long-term applications.

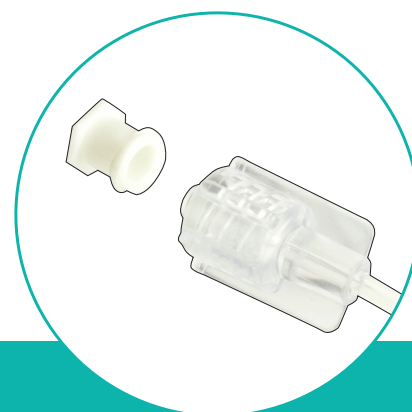

Available gas sampling lines



Nasal Sampling Line

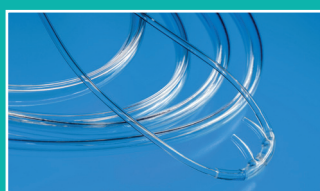

Nasal Cannula



Nasal/Oral Cannula




Airway Adapter Sampling Line

Reduction of contamination risks







The integrated bacteria filter of the VersaStream® gas sampling lines prevent contamination of gas-conducting components in and at the patient monitor. Transmission of germs, e.g. by contamination at the sample gas connection port of the monitor or at the water separation unit is prevented.

Water traps or water absorbers that are used in many gas monitoring applications in anaesthesia or intensive care – require either no replacement; or reduced need for replacement and emptying when using VersaStream® gas sampling lines.

Alarms triggered at the gas monitor by full water traps are significantly reduced. The associated patient monitoring labour and consumables costs are reduced considerably.

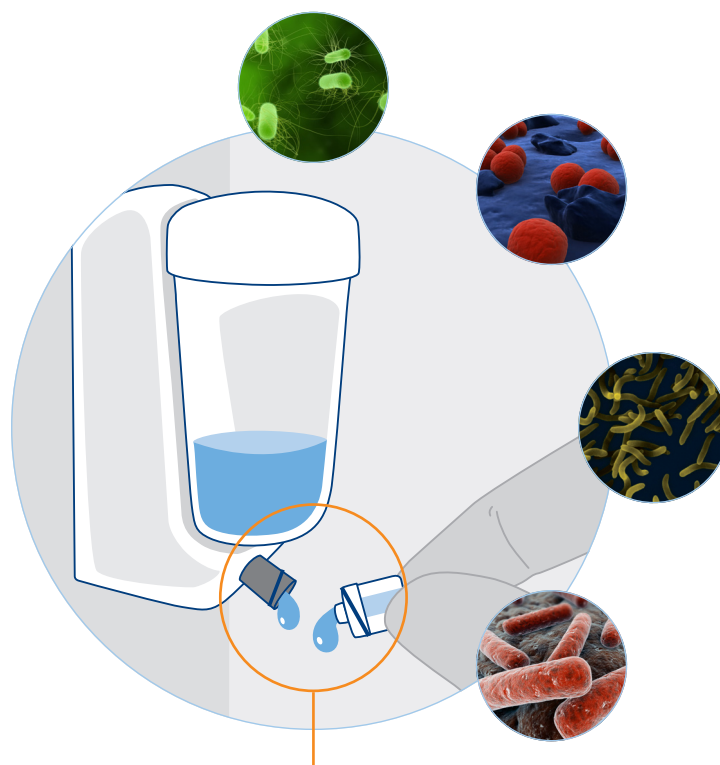

Contamination risk where standard lines are used with systems that employ water traps. VersaStream® lines resolve this potential infection control issue.

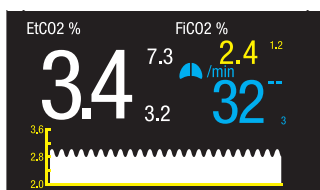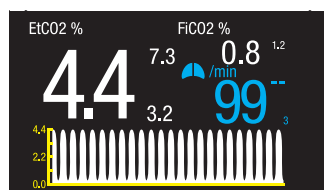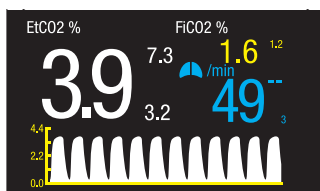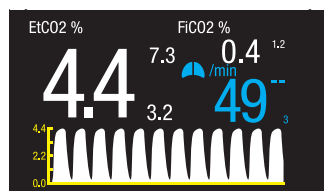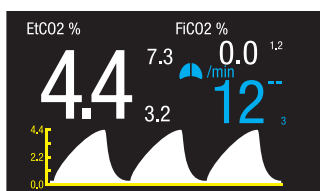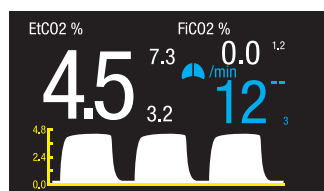

VersaStream®

LDS filter – Low-dead-space bacteria filter

Reliable measured values – even in cases of high breathing or ventilation frequency

The dead space volume of the filter unit has been reduced to only a few micro litres, thus facilitating precise measurement of CO₂ concentrations in the breathing gas – even at high respiration rates.



Diverse areas of application

VersaStream® gas sampling lines are suitable for use in areas, such as:

- Short and long term patient monitoring applications
- Monitoring of intubated and non-intubated patients
- Emergency medicine
- Patient transport
- Anaesthesia
- Intensive-care medicine – for adults and children
- Sedation
- Pain therapy, pain therapy management
- Homecare
- Sleep clinics and CPAP-ventilation



Long-term application – Nasal Cannula		Packaging unit	Viamed Part Number	Model Number
CO ₂	Adult	25 pcs./box	4420857	NA-L-LLM-200
	Paediatric	25 pcs./box	4420858	NP-L-LLM-200
	Infant	25 pcs./box	4420859	NI-L-LLM-200
O ₂ /CO ₂	Adult	10 pcs./box	4420860	NAO2-L-LLM-200
	Paediatric	10 pcs./box	4420861	NPO2-L-LLM-200
	Infant	10 pcs./box	4420862	NIO2-L-LLM-200

Long-term application – Nasal/Oral Cannula		Packaging unit	Viamed Part Number	Model Number
CO ₂	Adult	25 pcs./box	4420867	ONA-L-LLM-200
	Paediatric	25 pcs./box	4420868	ONP-L-LLM-200
O ₂ /CO ₂	Adult	10 pcs./box	4420869	ONAO2-L-LLM-200
	Paediatric	10 pcs./box	4420870	ONPO2-L-LLM-200

Long-term application – Airway Adapter Sampling Line		Packaging unit	Viamed Part Number	Model Number
Adult/Paediatric		25 pcs./box	4420873	A22-L-LLM-200
Infant/Neonatal ID ≤ 4.0 mm		25 pcs./box	4420874	A4.0-L-LLM-200

Long-term application - Sampling Line		Packagin unit	Viamed Part Number	Model Number
Luer Lock sampling line - LLM to LLM		25 pcs./box	4420876	LLM-L-LLM-200
Airway Adapter only with Luer Lock (LLF) connection If required, for use with 4420876		25 pcs./box	4420531	A22-LLF

Sampling flow rate ≤ 230ml/min

Long term ≤ 72 hrs

LLM Luer Lock Male – general monitoring applications that have a female Luer Lock connection port

LLF Luer Lock Female for gas sampling lines with LLM connector

O₂/CO₂ = O₂ delivery and CO₂ monitoring
CO₂ = CO₂ monitoring only

CE/CE0482

Specifications subject to change



Distributed by:



Viamed Limited · 15 Station Road · Cross Hills
Keighley · West Yorkshire · BD20 7DT · United Kingdom
Tel: +44 (0)1535 634 542 Fax: +44 (0)1535 635 582
Email: info@viamed.co.uk Website: www.viamed.co.uk



Part Number: 4490101
Date: 10/18

