


UNIVERSAL MOUNTING CLAMP INSTRUCTIONS
P/N CP2345

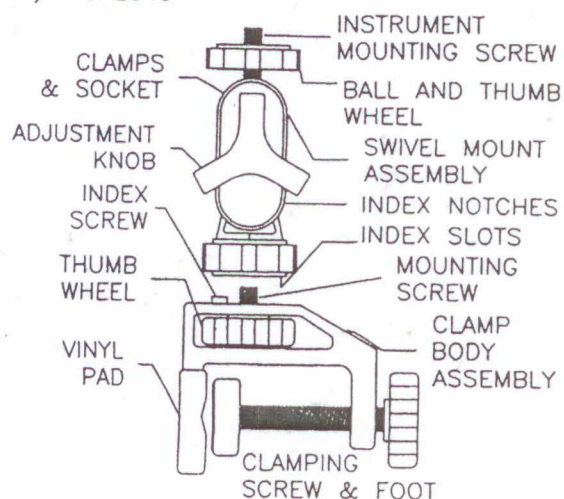

The Universal Mounting Clamp is used to attach the AX300 Oxygen Analyzer or MX300 Oxygen Monitor to IV poles, ventilators and other objects up to 1 ½ inch thick. The Mounting Clamp offers 360-degree rotation and 180-degree movement in all directions providing maximum flexibility in selecting a mounting site.

The clamp is supplied in two parts. The Clamp Body Assembly and Swivel Mount Assembly. To attach the Universal Clamp to the Swivel Mount use the Mounting Screw and Thumb Wheel shown in the diagram above.


Slots in the base of the Swivel Mount allow 360-degree rotation in 45-degree increments. When the desired orientation is obtained, insert the index screw into the slot and tighten the thumb wheel.

To attach the instrument to the ball mount loosen the Adjustment Knob until the ball mount screw is loose. Insert the Instrument Mounting Screw into the brass nut located on back of the AX/MX300 and tighten using the Ball Thumb Wheel.

Mount the Universal Clamp in the desired location and tighten the Clamping Screw. Adjust instrument to desired position and tighten the Adjustment Knob to secure the instrument in place.

A78141 Rev.0

UNIVERSAL MOUNTING CLAMP INSTRUCTIONS
P/N CP2345

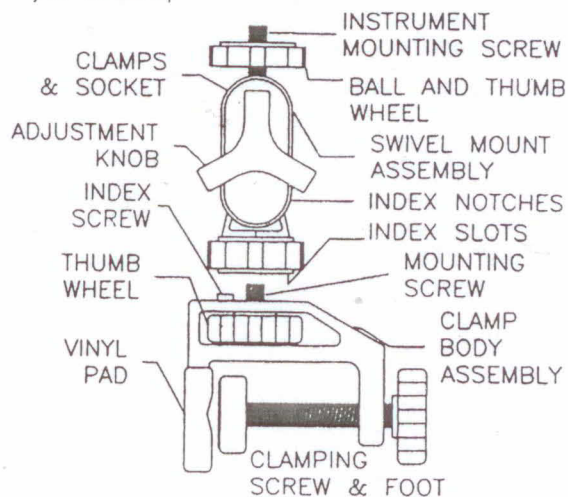

The Universal Mounting Clamp is used to attach the AX300 Oxygen Analyzer or MX300 Oxygen Monitor to IV poles, ventilators and other objects up to 1 ½ inch thick. The Mounting Clamp offers 360-degree rotation and 180-degree movement in all directions providing maximum flexibility in selecting a mounting site.

The clamp is supplied in two parts. The Clamp Body Assembly and Swivel Mount Assembly. To attach the Universal Clamp to the Swivel Mount use the Mounting Screw and Thumb Wheel shown in the diagram above.

Slots in the base of the Swivel Mount allow 360-degree rotation in 45-degree increments. When the desired orientation is obtained, insert the index screw into the slot and tighten the thumb wheel.

To attach the instrument to the ball mount loosen the Adjustment Knob until the ball mount screw is loose. Insert the Instrument Mounting Screw into the brass nut located on back of the AX/MX300 and tighten using the Ball Thumb Wheel.

Mount the Universal Clamp in the desired location and tighten the Clamping Screw. Adjust instrument to desired position and tighten the Adjustment Knob to secure the instrument in place.

A78141 Rev.0

P/N - 0121200