




Tom Thumb Infant Resuscitator

This easy to use, compact and robust infant resuscitator, has been designed to enable safe, controlled and fatigue-free resuscitation

- Compact and easy to use
- Complete control over PIP and PEEP
- Safe and consistent resuscitation
- Fatigue-free operation

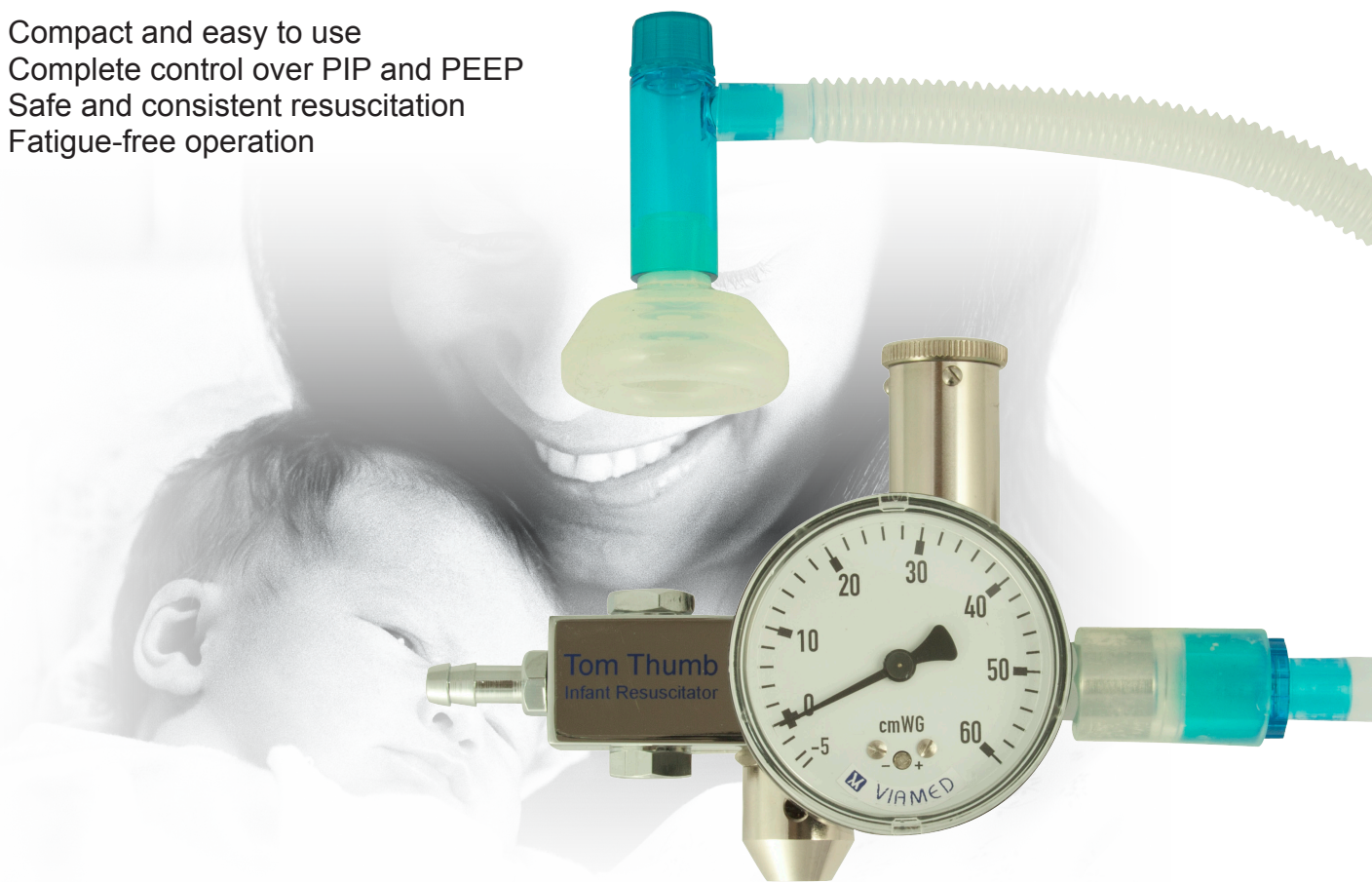

The Tom Thumb Infant Resuscitator has been designed to replace manual ventilation by the bag and mask method, which can often become tiring and inconsistent. The resuscitation is carried out by using the T-piece occlusion method, thus providing complete control and increased safety over infant resuscitation.

The NeoPEEP neonatal resuscitation circuit, supplied with the Tom Thumb, incorporates a variable PEEP valve and facilitates the delivery of safe and consistent breathing patterns to neonates; by enabling the delivery of consistent PEEP to assist in maintaining functional residual capacity.

Tom Thumb Infant Resuscitator

The Tom Thumb is available in several configurations

Part Number

Configuration One: TT480 for connection to an existing flowmeter or air/oxygen blender

0310030

Includes:

Variable blow off valve

Fixed blow off valve (predefined, e.g. at 45cm H₂O)


Integral medirail mounting clamp. Pole mounting clamp also available

NeoPEEP; variable PEEP neonatal resuscitation circuit (single patient use)

Single use silicone round facemask, size 1

Note: This version of the Tom Thumb has a 'Christmas Tree' inlet for the connection to an external flowmeter/blender via oxygen tubing

Specification subject to change

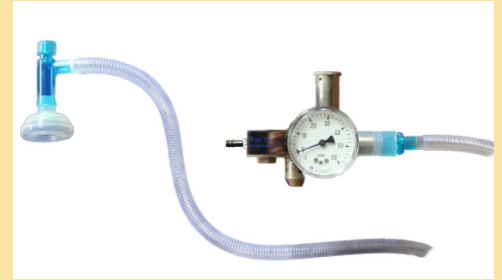

TT480 Infant Resuscitator



Rear of Tom Thumb with rail clamp

Replacement Parts and Accessories

NeoPEEP neonatal resuscitation circuit (single patient use)

3210011

Adapter: 15 mm I.D. - 15mm I.D.

0120140

Reusable; for connecting the NeoPEEP to the Tom Thumb

Adapter: 15 mm I.D. - 15 mm I.D. *Pack of 20*

0120151

Disposable; for connecting the NeoPEEP to the Tom Thumb

Single use silicone round facemask - size 0

3210070

40 mm (Height) x 27.4 mm (Inner Width) x 48.6 mm (Diameter) - Connector Diameter 15 mm O.D.

Single use silicone round facemask - size 0.5

3210067

44.4 mm (Height) x 32 mm (Inner Width) x 54 mm (Diameter) - Connector Diameter 15 mm O.D.

Single use silicone round facemask - size 1

3210071

47.9 mm (Height) x 40.5 mm (Inner Width) x 60 mm (Diameter) - Connector Diameter 15 mm O.D.

0380000

Annual service - Tom Thumb

Service and calibration of instrument to original specification.

Includes all 'O' rings, additional parts are chargeable.

Calibration certificate included.

Other products in our Infant Resuscitation range include:

Infant Resuscitation Cabinet

Apgar Timer

Air/Oxygen Blenders

We also supply a range of equipment for use during phototherapy, infant warming and oxygen monitoring.

CE 0086

Distributed by:



Viamed Limited · 15 Station Road · Cross Hills
Keighley · West Yorkshire · BD20 7DT · United Kingdom
Tel: +44 (0)1535 634 542 Fax: +44 (0)1535 635 582
Email: info@viamed.co.uk Website: www.viamed.co.uk



Part Number: 0390004
Date: 11/17