

HOSPITALS

Displays measurement of  
breath alcohol of a patient

Department Manager  
«DEPARTMENT»  
«COMPANY»  
«ADDRESS1»  
«ADDRESS2»  
«ADDRESS3»  
«CITY» «STATE»  
«ZIP»



Dear Sir or Madam,


**AlcoTrue M Breathalyser provides rapid and accurate  
determination of breath alcohol concentration**

The AlcoTrue M is a professional handheld medical breathalyser<sup>intended</sup> for quick and accurate determination of breath alcohol concentration in the exhaled breath of patients with clinical suspicion of alcohol intoxication, ~~which rivals its competitors and closely matches evidential devices~~.

~~The AlcoTrue M is designed for diagnostic use e.g. hospitals, rehabilitation clinics and other medical scenarios. They are used when there are concerns regarding blood alcohol level or alcohol consumption, for example in: emergency medicine, psychiatry, addiction clinics or at a methadone maintenance clinic to avoid side effects.~~

The main benefits include:

- BOLD TYPE** — {
- Fast warm-up time < 6 seconds; allowing for quick and easy diagnosis.
  - Rapid, precise and reliable clinically accurate measurements ~~even at sub-zero ambient conditions, can be used on site or in the field.~~
  - Easy to read display.
  - This high spec product provides up to five measurement modes, ideal for any situation;
    - **Active mode** – To determine the precise alcohol level in the end-expiratory breath of the person to be tested after achieving a minimal breath volume from 1.2 L. It is ensured that only deep lung breath (alveolar breath) is analysed.
    - **Fast mode** – Intended for persons suffering from COPD, Asthma or reduced lung capacity, whom might not be able to provide the standard minimum sample volume.
    - **Passive mode\*** – To obtain a precise measurement, the sample uptake must be triggered when the person is breathing out. The measurement mode is to determine the quantitative value of alcohol concentration.
    - **Scan mode\*** – Qualitative testing of the ambient air alcohol content, qualitative determination of alcohol content in exhaled air of a person who is not able to provide a breath sample (e.g. unconscious person), qualitative test of "non-alcoholic" drinks for alcohol content.
    - **Mouth alcohol mode\*** – Consists of two tests with a determined waiting period between first and second measurement.

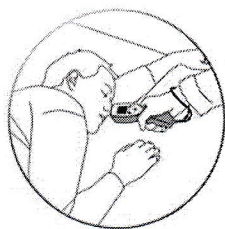

\* Measurement modes only available if activated.

- Also available with a portable printing system.

Please find enclosed further information regarding the AlcoTrue M. Should you wish to discuss this in further detail please do not hesitate to contact a member of our sales team on 01535 634542 or by emailing [alcotrue@viamed.co.uk](mailto:alcotrue@viamed.co.uk).

Yours faithfully

Catrin

Catrin Hollings  
Marketing

Should you no longer wish to receive  
marketing material via post,  
sign-up to our email mailing list at  
[www.viamed.co.uk/sign-up](http://www.viamed.co.uk/sign-up)