



Derek Lamb <liquidgands@gmail.com>

FWD: WG: AW: oxytrue as housing

3 messages

Steve Nixon <steve.nixon@viamed.co.uk>
To: Derek I Lamb <derek.lamb@viamed.co.uk>

24 July 2013 13:16

Drawing attached

----- Forwarded Message -----

From: "Lindner, Bernd" <b.lindner@bluepoint-medical.com> (By way of info@viamed.co.uk)
To: SNIXON (Steve Nixon) <GoldMine User>
Date: Wed, 24 Jul 2013 13:16:13 +0100
Subject: WG: AW: oxytrue as housing

Hello Steve,

I think shape 1 looks best. Can you agree. We want to push the project forward.
With kind regards
Bernd

Von: Radev, N.
Gesendet: Dienstag, 23. Juli 2013 09:35
An: william.pan@justec.cn
Cc: Lindner, Bernd; Froehlich, Heidi
Betreff: AW: AW: oxytrue as housing

Hi William,.

We believe the optimal location and shape for this extra warning LED is a shape 2 as attached picture.

It looks good and it will be easy for implementation. If there is no place on the main pcb then this led can be mounted on an extra small pcb and fit on the housing (melted pins) from inside.

Best regards

Ned

Dr.-Ing. Nedyalko Radev

/Development Engineer/

bluepoint MEDICAL GmbH & Co. KG

Tel +49-(0)-38823-548-8024
Fax +49-(0)-38823-548-8029
Mail n.radev@bluepoint-medical.com
Web www.bluepoint-medical.com

Trade register: Schw erin, Reg.-No.: HRA 2600

Disclaimer: This message may contain confidential and /or privileged information. If you are not the addressee or authorized to receive this for the addressee, you must not use, copy, disclose or take any action based on this message or any information herein. If you have received this message in error, please advice the sender immediately by reply e-mail and delete this message.

Thank you for your cooperation.

Von: williampan [mailto:william.pan@justec.cn]

Gesendet: Dienstag, 23. Juli 2013 05:41

An: Radev, N.; Froehlich, Heidi; Lindner, Bernd

Betreff: Re: AW: oxytrue as housing

Hello Ned,

How about this design?

ãã

Best Regards,

æ½«æ±

ä½3;æå¼·ç§æ(æ±3;)æé,
æ±3;å,åå±å°©®¶e¼åæ°·88å·éé¼å·¥ä,å°64æ 2æ¼¼

William Pan

Justec Shenzhen Co., Ltd

2F,Building 64,Majialong Industry Zone, Shenzhen,China,518054

TEL:+86-755ï¼26068931

Fax:+86-755-26068933

Mobile:+86-13631608890

www.justec.cn

2013-07-23

===== 2013-07-22 16:20:39 ææ¥ä¼ä,åé¼¼ =====

Hi William,

Water proof is not a problem.

You should set design to have transparent part (slightly milky colour) as cover for this additional warning LED. This part has to be from ABS material and will be glued to the Housing with solvent glue or UV glue.

Using the keypad back light will be not sufficient enough on daylight.

Best regards

Ned

Dr.-Ing. Nedyalko Radev

/Development Engineer/

bluepoint MEDICAL GmbH & Co. KG

Tel +49-(0)-38823-548-8024

Fax +49-(0)-38823-548-8029

Mail n.radev@bluepoint-medical.com

Web www.bluepoint-medical.com

Trade register: Schw erin, Reg.-No.: HRA 2600

Disclaimer: This message may contain confidential and /or privileged information. If you are not the addressee or authorized to receive this for the addressee, you must not use, copy, disclose or take any action based on this message or any information

Updated drawings attached

----- Forwarded Message -----

From: "Lindner, Bernd" <b.lindner@bluepoint-medical.com> (By way of info@viamed.co.uk)
To: SNIXON (Steve Nixon) <GoldMine User>
Date: Wed, 24 Jul 2013 13:17:23 +0100
Subject: WG: AW: AW: oxytrue as housing

Here the shape 1 version.

Kind regards
Bernd

Von: williampan [william.pan@justec.cn]
Gesendet: Mittwoch, 24. Juli 2013 10:50
An: Lindner, Bernd; Radev, N.
Cc: Froehlich, Heidi
Betreff: Re: AW: AW: oxytrue as housing

Hello Bernd,

There is no space for the shape 2.
Here is the design based on shape 1.
ãã

Best Regards,

æ½å«æ±

ä½▷3;æå¾·ç§æ(æ·±å▷3;)æéå¬å,
æ·±å▷3;å,åå±±±å°é©¬å®¶¶e¾å¤§æ°è·88å·éé¾å·¥ä,å°64æ 2æ¥¾

William Pan
Justec Shenzhen Co., Ltd
2F,Building 64,Majialong Industry Zone, Shenzhen,China,518054
TEL:+86-755ï¼26068931
Fax:+86-755-26068933
Mobile:+86-13631608890
www.justec.cn
2013-07-24

===== 2013-07-24 14:44:00 æ¬äæ¥ä¿,åéï¼ =====

Hello William,

shape 1 looks best.

With kind regards

Bernd

[Quoted text hidden]

[Quoted text hidden]

This message has been scanned for malware.

4 attachments




ap_mobile.jpg
39K



ap_mobile-2.jpg
44K



ap_mobile-3.jpg
42K

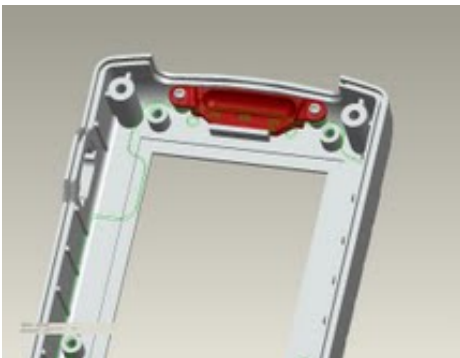

ap_mobile-4.jpg
68K

Steve Nixon <steve.nixon@viamed.co.uk>
To: Derek I Lamb <derek.lamb@viamed.co.uk>

24 July 2013 13:18

----- Forwarded Message -----

From: "Lindner, Bernd" <b.lindner@bluepoint-medical.com> (By way of info@viamed.co.uk)
To: SNIXON (Steve Nixon) <GoldMine User>
Date: Wed, 24 Jul 2013 13:18:08 +0100
Subject: AW: AW: oxytrue as housing

We only want to add one additional alarm light, which is shape 1.

With kind regards
Bernd

Von: steve.nixon.viamed@googlemail.com [steve.nixon.viamed@googlemail.com] im Auftrag von Steve Nixon
[steve.nixon@viamed.co.uk]

Gesendet: Mittwoch, 24. Juli 2013 10:52

An: Lindner, Bernd

Betreff: Re: AW: oxytrue as housing

Bernd

Why do we need the extra alarm indicator/s? Is it legislation, because the TFT is not as clear as OLED, or just to provide extra all around indication?

Are you planning to fit 3 warning lights (one on top and one either side) as per the drawings or just one on top?

If 3 lights â then I suggest shape 2 for the top.

If 1 light â then perhaps shape one is the better option.

Steve

On 24 July 2013 07:43, Lindner, Bernd <

b.lindner@bluepoint-medical.com> wrote:

Hello Steve,

I think shape 1 looks best. Can you agree. We want to push the project forward.

With kind regards

Bernd

[Quoted text hidden]

[Quoted text hidden]

[Quoted text hidden]

www.justec.cn

[Quoted text hidden]

[Quoted text hidden]

—

Steve

Steve Nixon
Director - Viamed Ltd.

This message has been scanned for malware.