


## Ceratherm 600-2/3 Radiant Warmers

### Installation Instructions for modified casing (part number 521-7198-3)

The new casing (part number: 521-7198-3) has six holes, each with a diameter of 8 mm (previously 6 mm).

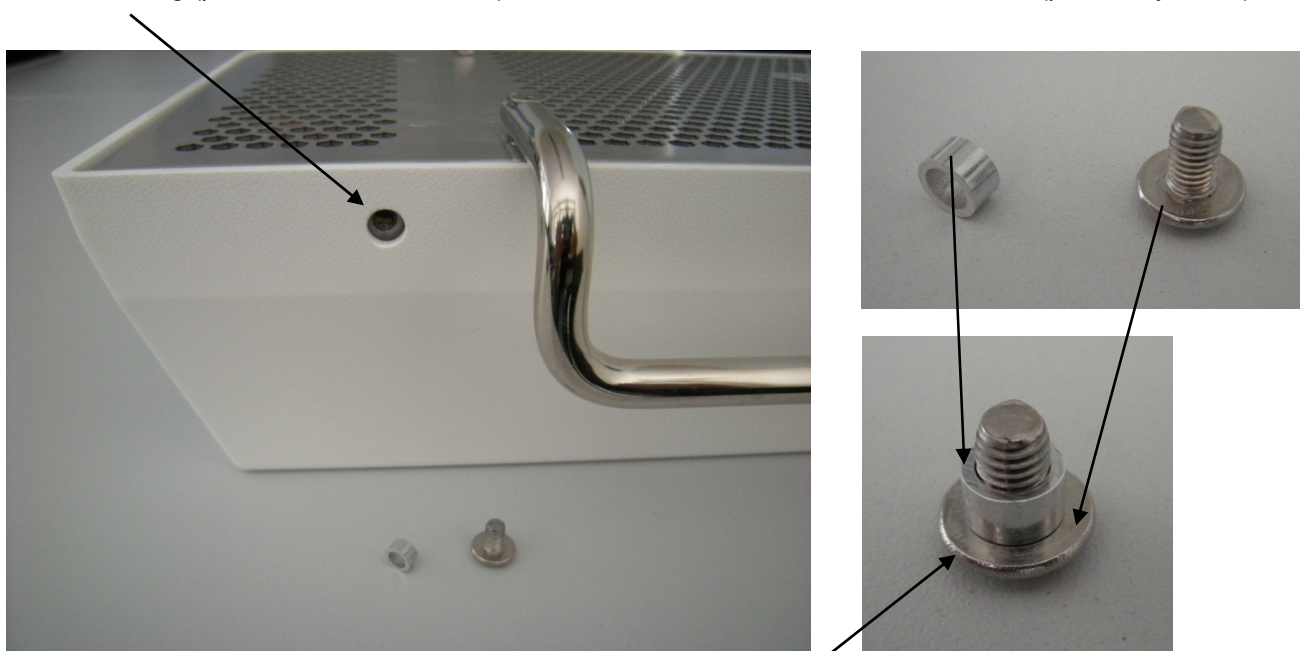

The existing six mounting screws (M5 x 8) should be assembled with the supplied 7.0 x 1.0 x 4.0 support sleeves (part number: 521-79940100.10, set of 6).

Each screw should then be fitted, together with a support sleeve.



Using a screwdriver, tighten by hand until the screw stops. The screw head should be tightened against the supporting sleeve, not against the casing.



2

### Caution!

**Ceratherm 600-2:** Affix the new minimum distance warning labels (part number 521-GF2704.07N, pack of two labels) and reaffix the existing type plate removed from the old casing. Use the old casing as a guide for the label placement.

**Ceratherm 600-3:** Affix a minimum distance warning label (part number 521-GF2704.07N, pack of 2 labels) to each side of the casing as shown below. NOTE: The type plate is located on the rear of warmer chassis.



Dispose of the old casing in accordance with local regulations.

\* \* \*