

breathalyser



AlcoTrue[®] P Breathalyser

Professional breathalyser for rapid accurate determination of breath alcohol concentration



Ideal for use within

- Police
- General Workplaces
- Prisons & Customs
- Transport: Road/Rail/Air



Professional Breathalyser for rapid accurate determination of alcohol concentration in exhaled air.



AlcoTrue P Standard Package

- AlcoTrue P Breathalyser; with active & fast modes**
- 25 mouthpieces (individually wrapped)
- Battery holder (4 x AA)
- 4 x AA batteries
- Instructions for use
- Small hard carrying case

** Further modes can be activated by Viamed, please confirm when ordering.

Key Features of the AlcoTrue P Breathalyser

- Professional fuel cell technology
- High precision intelligent microprocessor sampling system
- Automatic self-test after switch on
- Rapid, precise and reliable measurements, even at sub-zero ambient conditions
- Intuitive user interface
- Ready to test < 6 seconds, even at temperatures below zero
- Results are shown on a 1.77" colour display
- Compact, ergonomic design
- Low maintenance and spare part costs

Measurement Modes

The **AlcoTrue P** provides up to five measurement modes, ideal for any situation:

■ Active mode AM

To determine the precise alcohol level in the end-expiratory breath of the person to be tested after achieving a minimal breath volume from 1.2 L. It is ensured that only deep lung breath (alveolar breath) is analysed.

■ Fast mode FM

Intended for persons suffering from COPD, Asthma or reduced lung capacity, whom might not be able to provide the standard minimum sample volume.

■ Passive mode PM*

To obtain a precise measurement, the sample uptake must be triggered when the person is breathing out. The measurement mode is to determine the quantitative value of alcohol concentration.

■ Scan mode SM*

Qualitative testing of the ambient air alcohol content, qualitative determination of alcohol content in exhaled air of a person who is not able to provide a breath sample (e.g. unconscious person), qualitative test of "non-alcoholic" drinks for alcohol content.

■ Mouth alcohol mode MA*

Consists of two tests with a determined waiting period between first and second measurement.

* Measurement modes only available if activated.



Applications

The AlcoTrue P Professional Breathalyser is ideal for use within:

- Police
- General Workplaces
- Prisons & Customs
- Transport: Road/Rail/Air

AlcoTrue P Breathalyser Kit

- AlcoTrue P Breathalyser
- Battery holder (4 x AA)
- 4 x AA batteries
- Li-Poly battery
- UK or EU* main power adapter with USB data cable
- Vehicle power adapter
- 25 hygienic mouthpieces (individually wrapped)
- AlcoPrint P printer
- Printer cable
- AlcoPrint UK or EU* mains power cable
- AlcoPrint vehicle power cable
- Single roll of printer paper
- PC software & USB cable
- Large hard carry case

* Must be specified at time of ordering.



The AlcoPrint P Printer

The **AlcoTrue P** Breathalyser Kit includes a AlcoPrint P USB wireless printer, this can print wirelessly from a range of 5 - 10 metres or connected directly to the unit with the USB cable provided. The AlcoPrint P printer is powered by either a vehicle or mains* power cable (provided). The AlcoPrint P printer is a thermal transfer printer.

The Printer Set is also available to purchase separately and includes: AlcoPrint P Printer, Printer Cable, UK or EU* Mains and Vehicle adapter and a roll of printer paper.




AlcoPrint USB Wireless Printer

PC Software Kit

The **AlcoTrue P** PC Software Kit includes a USB cable and PC software CD.

The PC Software CD is also available without the USB cable.

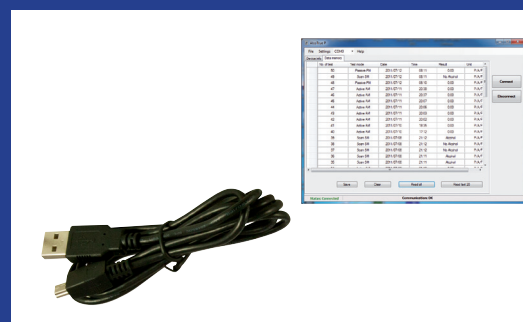

PC Software

Additional Accessories

Mouthpieces are available in packs of 25.

AlcoPrint paper rolls are available in packs of 5.



Clinically tested in-built non-return valve and water trap


AlcoTrue® P Breathalyser

CE0482


Sensor type		Professional Electrochemical sensor	
Sampling system		High precision sampling system 0.2ml	
Measurement modes		Active mode (default) Fast mode (activated) Passive mode (activated) Mouth alcohol mode (if activated) Scanning mode (if activated)	
Measurement range			
Breath alcohol concentration (BrAC)		0.00 to 2.6 mg/l, 0 to 2600 µl, 0.000 to 0.550 g/210L 0 to 260 µg/100mL, 0 to 260 µg% (Ethanol mass per breath volume at 34°C and 1013 hPa.)	
Blood alcohol concentration (BAC)		0.00 to 5.50‰, 0.000 to 0.550% , 0.00 to 5.50 g/l 0 to 550 mg/100ml, 0 to 550 mg%, 0.00 to 5.50 mg/ml (Ethanol mass per blood volume resp. mass at 20°C and 1013 hPa)	
This data is given using a conversion factor of 2100:1 between breath - and blood alcohol concentrations. Other conversion factors are available upon request.			
Ambient conditions			
For operations		-5 to 50°C (23 to 122°F), 10 to 100 % R.H. (non-condensing), 600 to 1400 hPa	
For storage		-20 to 65°C (-4 to 149°F), 10 to 100 % R.H. (non-condensing), 600 to 1400 hPa	
Sample conditions			
Minimal volume		Adjustable, 1.2L (set as default)	
Measurement accuracy			
Reproducibility with ethanol standard			
0 to 0.476 mg/l	0.008 mg/l	0 to 1‰	0.017 mg/l
> 0.476 mg/l to 1 mg/l	4% of measured value	> 1‰ to 2‰	4% of measured value
> 1 mg/l	8% of measured value	> 2‰	8% of measured value
Sensitivity drift of measured value/month		Typically 0.5%	
Ready for operation: After switching on		< 6 seconds	
Configurable Units		µg/l µg/100ml % g/100ml mg% P/A/F	g/210l µg% g/l mg/100ml mg/ml P/F
Dimensions (W x H x D)		185 mm x 55 mm x 40 mm	
Weight		240g (including Alkaline batteries) 195g (including Li-Poly battery)	
Display		1.77" TFT Colour display	
Data memory		Up to 35,000 test result measurements	
Standards		EN 15964:2011 certified, EN 55011:2007, EN 61326:2006, IEC 60068-2-6:2007, IEC 60068-2-27:2008, IEC 60068-2-64:2008, IEC 68-2-32:1975 + A1/A2:1980, EN 60529:1991 +A1:2000	
Calibration required		6 months	
Part Numbers		AlcoTrue P Standard Package AlcoTrue P UK Kit AlcoTrue P EU Kit	4710011 4710032 4710033

Screen Displays

Qualitative results indications (Scanning and Mouth Alcohol modes)

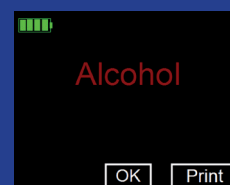


Alcohol detected in gas sample

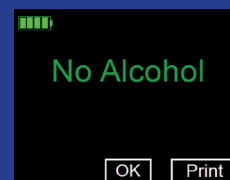


No alcohol detected in gas sample


Quantitative results indications (Active, Fast, Passive and Mouth alcohol modes)

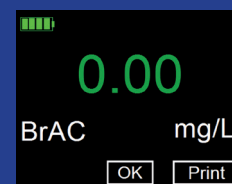

Pass

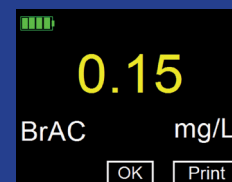

Alert

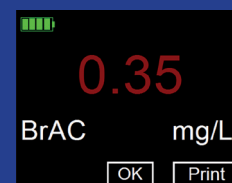

Fail

Specification Subject to Change

Distributed by:



Angel Healthcare Solutions Ltd
Ponderosa Business Park
Smithies Lane
Heckmondwike, West Yorkshire
WF16 0PN

Tel - +44 (0)1924 412133
Fax - +44 (0)1924 414829
info@angelhealthcaresolutions.com



Viamed Limited · 15 Station Road · Cross Hills
Keighley · West Yorkshire · BD20 7DT · United Kingdom
Tel: +44 (0)1535 634 542 Fax: +44 (0)1535 635 582
Email: info@viamed.co.uk Website: www.viamed.co.uk



Part Number: 4790022
Date: 11/16